

# LESCOLTON – POWERED BY ELOS ME

## USER GUIDE ICECOOL I DUAL SYSTEM



390 Market St 200  
San Francisco, CA, 94102

Questions? Get in touch with us by live chat or email us:  
[info@lescoltonglobal.com](mailto:info@lescoltonglobal.com)

## CONTENTS

Product Introduction .....	03
Contraindications.....	03
Considerations.....	03
Safety warning.....	03
Safety operation regulations .....	03
Product structure .....	04
How to use the device.....	05
Operation suggestions.....	06
Treatment intervals .....	06
Specifications .....	06
Table of simple troubleshooting .....	07
Warranty registration .....	07

## **INTRODUCTION**

- The new Lescolton powered by Elos Me is a dual system with the functions of hair removal and skin rejuvenation. It integrates the versatile application of optical e-light technology in the field of Beauty & Health.
- Before using the device, please read the instructions carefully.
- Suitable for both men and women.
- Individuals with dark skin should start with the lowest level and then adjust to the proper level based on their skin's adaptability.

## **CONTRAINDICATIONS**

- Not suitable for individuals under 16 years old.
- Avoid use in cases of skin inflammation.
- Avoid use in cases of photosensitivity.
- Avoid use if you have an implanted medical electronic device that may lead to operational errors.

## **CONSIDERATIONS**

- This product has not been tested on pregnant women.
- Avoid using the product directly on tattoos as they may absorb a significant amount of energy.
- Refrain from using around the eyes, neck, or nipples as the skin in these areas is relatively sensitive.
- If skin remains inflamed or swollen after 24 hours, discontinue use until symptoms subside.


## **SAFETY WARNING**

- Do not block the vents during use.
- Ensure that no water or moisture enters the device.
- Using the device in inflammable or explosive conditions is prohibited.

## **SAFETY OPERATION REGULATIONS**

- Avoid directing the light at the eyes.
- Do not look directly at the light.
- Before use, ensure that you have removed all metallic objects such as rings and bracelets to prevent device malfunction.
- Operate the machine in a room with a temperature below 25 degrees Celsius.

# PRODUCT STRUCTURE

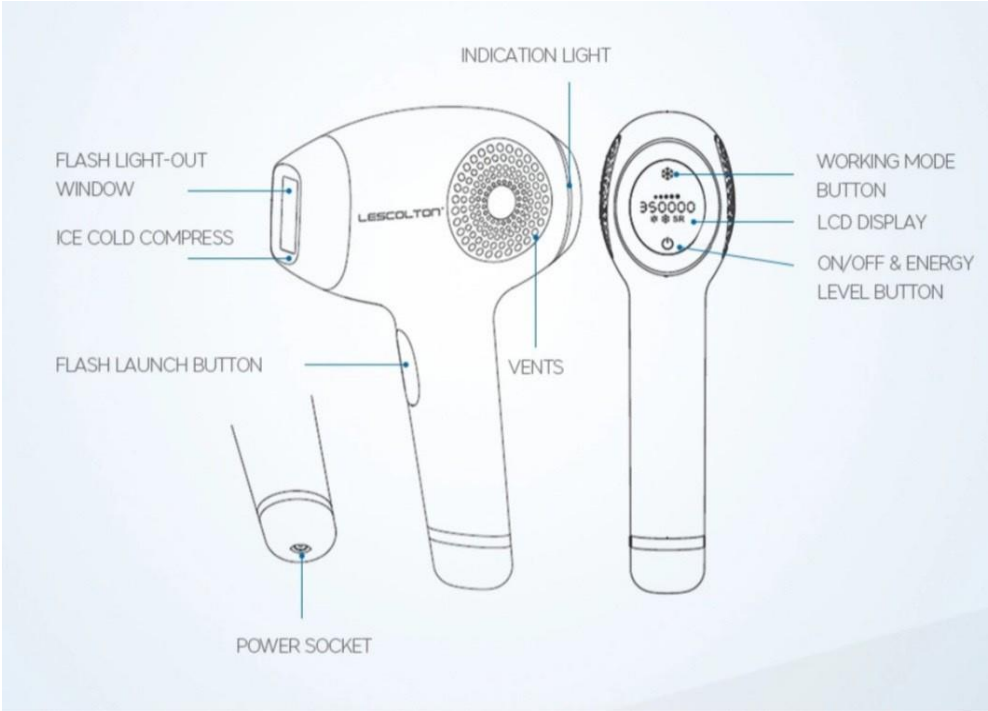

## HOW TO USE

- Preferably, wax all the hair around the treatment area before use. If you can't wax, you can shave the treatment area.
- Insert the power cable into the machine and plug it in.
- Turn on the machine by pressing and holding the touch power button on the LCD screen for two or three seconds until the device turns on.
- Use the same power button to select the desired power level. Each press of the button will increase the power level. If you don't feel any tingling or discomfort at the lowest power levels, it's best to select the highest possible power level.
- Place the lamp on the skin, and once the green lights around the LCD screen of the device flash, press the fire button on the device handle. Repeat the shots on the same area 2 or 3 times for the best results. To trigger the next shot, wait for the green lights on the side of the machine to flash.
- For skin rejuvenation treatment, activate the SR mode. Once this mode is activated, simply place the lamp on the skin and press the fire button. This treatment mode can be used on the same body areas as the hair removal treatment.

### Treatment Modes

- Four treatment modes can be selected using the top button on the LCD screen: hair removal, hair removal with cooling, skin rejuvenation, and skin rejuvenation with cooling.
- Each time you press the top button, the treatment mode will change. When the top button lights up in blue upon pressing it, it means that the cooling mode is activated. When pressing the button displays the abbreviation "SR" on the LCD screen, it signifies that the skin rejuvenation treatment mode is active. If the "SR" disappears from the display upon pressing the button, it indicates that the hair removal mode is active.

### Automatic Use Mode

- To allow the machine to emit light pulses automatically without pressing the shutter button, switch to automatic mode. To do this, press the shutter button for four or five seconds until you hear two beeps, and the green lights around the LCD screen of the device will flash. The machine will automatically emit light pulses when the lamp is in contact with the skin. It is recommended to keep the lamp in contact with the skin until two or three pulses are emitted for faster results. To return to manual pulse emission mode, press the shutter button again for four or five seconds until you hear two beeps, and the green lights around the LCD screen of the device will flash.

## TREATMENT INTERVALS

Treat every three to five days during the first month, weekly during the second month, every 15 days during the third month, and every 30 days from the fourth month until permanent hair removal is achieved.

Permanent hair removal is a process that takes several months. Immediate results are not expected for two reasons:

1. It is necessary to wait until all body hair is in the active growth stage to be affected by the treatment, which may take several months.
2. Multiple treatments are required for hair in the active growth stage to weaken the hair root until it is permanently affected, preventing regrowth. Visible changes will occur over the course of several months.

Individuals with blond, grey, or red hair will typically require twice as many sessions as those with black and brown hair. Therefore, the treatment duration may extend to 6 to 12 months."

## **SPECIFICATIONS**

- Main Unit Size: 182x78x151mm
- Package (Box) Size: 231x191x93mm
- Power Supply: Input: 100-240V AC-50/60Hz
- Weight: 1.1 kg
- Irradiation Area: 4 cm<sup>2</sup>
- External Environment: Temperature Range: -20 to 50 degrees Celsius
- Warranty: 6-year warranty

## **TABLE OF SIMPLE TROUBLESHOOTING**

- No Power: Ensure the power cord is securely connected.
- Can't Turn On the Machine: Press the "on" button for 2 seconds.
- Level Indicator Is Flashing After Boot: The cartridge may not be properly installed.
- Automatic Shutdown: The machine shuts down automatically after 5 minutes of inactivity.