

ELOS ME & LESCOLTON

USER MANUAL INFINITY DUAL SYSTEM



390 Market St 200
San Francisco, CA, 94102

Questions? Get in touch with us by live chat or email us:
info@elosme-lescolton.com

CONTENTS

Product Introduction.....	03
Contraindications.....	03
Considerations.....	03
Safety warning.....	03
Safety operation regulations.....	03
Product structure.....	04
How to use it.....	05
Permanent hair removal treatment.....	05
Treatment intervals.....	05
Skin rejuvenation treatment.....	06
Specifications.....	06
Table of simple troubleshooting.....	06
Warranty registration.....	07

INTRODUCTION

- The new Lescolton powered by Elos Me is a dual system, with the function of hair removal and skin rejuvenation. It integrates the multiterm application of optical e-light technology in the field of Beauty & Health.
- Before use, please read the instructions carefully.
- Suitable for men and women.
- People with dark skin should start with the lowest level and then adjust to the proper level based on the adaptability of your skin.

Contraindications

- People under 16 years old.
- Inflammation of the skin.
- Photosensitivity.
- Implanted medical electronic device which may lead to operation errors.

Considerations

- The product has not been tested on pregnant women.
- Do not use the product directly on tattoos since they may absorb a large quantity of energy.
- Do not use around the eyes, neck or nipples since the skin of these areas is relatively vulnerable.
- If the skin continues to be inflamed or swollen after 24 hours, discontinue use until symptoms disappear.

Safety warning


- Do not block the vents during use.
- Do not allow the water or any moisture into the device.
- Use the device in inflammable or explosive conditions is prohibited.

Safety operation regulations

- Do not aim the light at eyes.
- Do not look directly at the light.

- Before use, please make sure that you have removed all metallic objects such as rings, bracelets and so on, to avoid malfunction of the device.
- Use the machine in a room with a temperature inferior to 25 degrees Celsius.

Product Structure

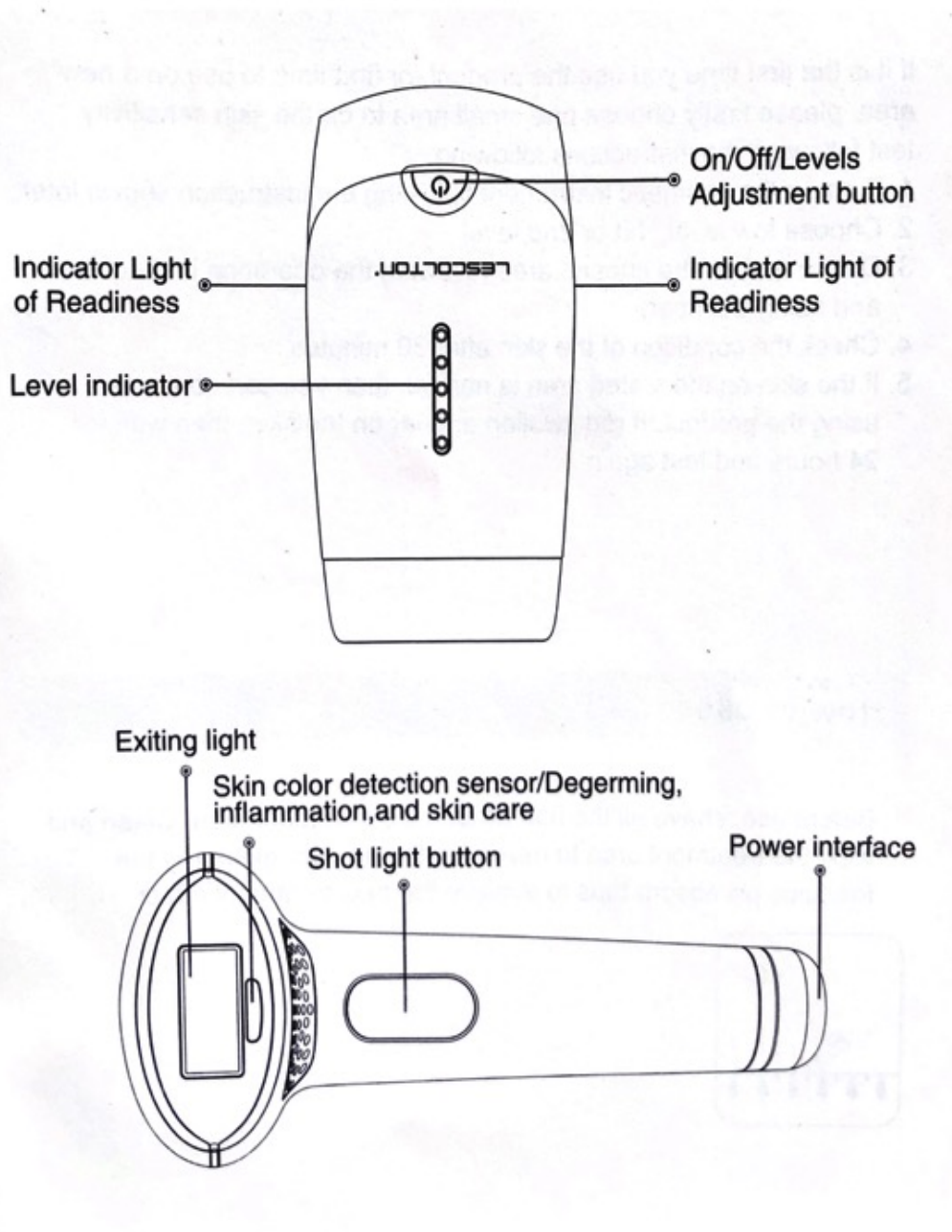

HOW TO USE

Permanent hair removal treatment

Before use, preferably wax all the hair around the treatment area. If you can't wax, you can shave the area.

Insert the power cable into the machine and plug it in.

Turn on the machine from the power by pressing it for 2 or 3 seconds until the device turns on.

With the same power button, select the desired power level. Each time the button is pressed, a power level will rise. If at the lowest power levels you do not feel tingling or some discomfort, it is best to select the highest possible power level.

Rest the lamp on the skin and once the green lights on the sides of the device flash, press the fire button on the device handle. Repeat shots on the same area 2 or 3 times for best results. To be able to issue the next shot it will be necessary to wait for the green lights on the side of the machine to flash.

To allow the machine to shoot automatically without pressing the shutter button, you can enter the automatic use mode. To enter automatic mode, press the shutter button for two or three seconds until two beeps are emitted. This way, the machine will emit the shots automatically when the lamp is resting on the skin. Keep the lamp supported until 2 or 3 shots are emitted for best results.

To return to manual shooting mode it will be necessary to press the shutter button again for 2 to 3 seconds until two beeps are heard again.

Treatment intervals

Treat every 3 to 5 days during the first month, weekly during the second month, every 15 days during the third month and every 30 days from the fourth month until you permanently remove the hair.

Definitive hair removal is a process of several months. The results are not seen immediately due to two reasons:

1- It is necessary to wait until all the body hair is in the active growth stage so that it can be affected by the treatment. For this, it is necessary to wait several months.

2- It is necessary to perform several treatments for hair that is in the active growth stage to weaken the hair root until it is permanently affected and the hair does not grow back.

The changes will be noticed over the months.

People with blond, grey or red hair will normally need twice as many sessions as people with black and brown hair, so the treatment can be extended to between 6 to 12 months.

Skin rejuvenation treatment

For this treatment it will be necessary to change the lamp.

With the machine turned off, pull the lamp to remove it and place the rejuvenation lamp.

Once the lamp is placed, follow the same steps to turn on the device, select the power level and shoot as explained previously for permanent hair removal.

Specifications

Size of the main unit: 182x78x151mm

Size of the package (box): 231x191x93mm

Power supply: input: 100-240V AC-50/60Hz

Weight: 1,1 kg

Irradation area: 4 cm²

External environment: Temperature: - 20 - 50 degrees Celsius

Warranty: 4 year warranty

Table of simple troubleshooting

- No power: The power cord is not well connected.
- Can't turn on the machine: Press the "on" button for 2 seconds.
- The level indicator is flashing after boot: The cartridge is not well installed,

- Automatic shutdown: The machine shuts down automatically after 5 minutes if it is not used.

Warranty registration

You can register your 4-year warranty period through this link:
<https://elosme.com/warranty.php>